

Alber Twion




Expand your possibilities and become more independent with the efficient Twion active drive. Compatible with virtually all common active wheelchairs, you can continue to use your current individually adapted wheelchair.

It only requires a small lightweight bracket to attach to your existing wheelchair and can be easily removed with a quick release mechanism when transporting or if you would prefer to reattach your regular manual wheels.

The patented EasyConnex charging plug magnetically connects to the charger socket (which also serves as an On/Off switch). The lithium-ion batteries in the wheel hubs guarantee years of trouble-free use and provide sufficient energy to handle going up and down steep slopes and allow you to cover greater distances.

- Motor, electronics and battery pack integrated in the wheel hub
- 24" wheels with Schwalbe Marathon Plus Evolution tyres
- Aluminium black powder-coated push rims
- Li-Ion batteries - allow up to 1,200 charge cycles and a range of approx 15km
- Battery charger - charge both batteries at same time - 4.5 hour charge time from flat
- Optional Mobility Plus App for Smartphone - additional features and modes

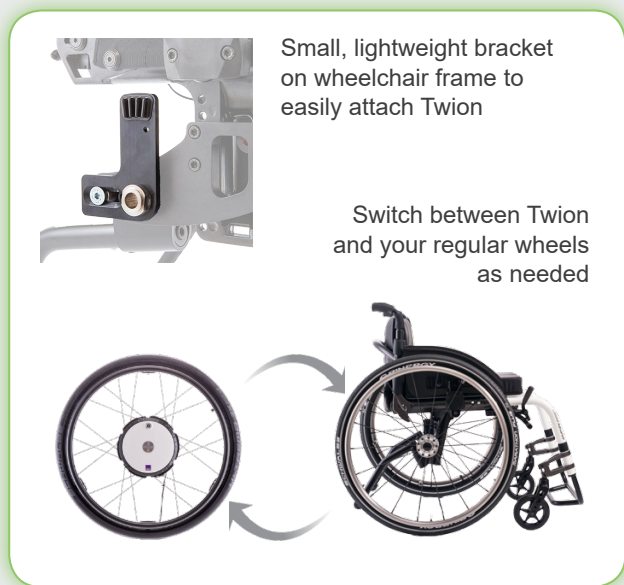

	Wheel Size	Stub Axle	Product Weight	Battery	Charge Time (from totally flat)	S.W.L.
Twion	24" (61cm)	12.3cm x 1.3 (stainless steel)	6.3kg (per wheel)	37V Li ion	4.5 hours	145kg